

1. DESCRIPTION

Equine herpesvirus type 1 (EHV-1) is a member of the Alphaherpesvirinae subfamily. EHV-1 is a major cause worldwide of epidemic abortion, perinatal mortality, respiratory disease and neurological disorders in horses. As with other herpesvirus infections, lifelong latency of the virus in the nervous system of the host occurs, and periodical reactivation of the virus can cause new outbreaks. Upper respiratory infection is the most common manifestation of EHV-1 infection. Commonly, young horses have depression, poor appetite, nasal discharge and cough. If young horses are housed or pastured together, many horses in the herd may experience disease from which they recover uneventfully. Disease may be mild or unapparent in older horses. Neurological signs occur infrequently as a result of EHV-1 infection; however, outbreaks of neurologic disease associated with fever, nasal discharge and cough have been reported in the USA and elsewhere. Neurologic signs may include incoordination that can progress to an inability to stand. Horses may be unable to urinate or may dribble small volumes of urine. Horses may have difficulty producing manure. Sometimes the neurologic signs are accompanied by cellulitis (inflammation or swelling of the limbs) and petechiae (small hemorrhages on the gums).

VetPCR™ EHV-1 Detection Kit is the direct detection of *Equine herpesvirus types-1* on the basis of a genetic database, so it can diagnose very fast and accurately. It can amplify only specific gene using the PCR (Polymerase Chain Reaction) method, and take only 3 hours for detection. Therefore, it is a very fast, accurate and reliable technique.

2. STORAGE

The components of VetPCR™ EHV-1 Detection Kit should be stored at -20°C. Under this condition, the kit is stable until expiration date stated on the label.

3. CONTENTS

	Kit 48	Kit 96	
VetPCR™ EHV-1 PCR Pre-mixture	48	96	tubes
DNase/RNase-free water	1	1	vial
EHV-1 PCR Positive control	1	1	vial
EHV-1 PCR Positive control Pre-mixture	4	8	tubes
Brig™ Molecular Weight marker	1	1	vial
Mineral Oil	1	1	vial
DNA purification kit (see step 6.1)	50	100	tests

4. SPECIMEN

Nasopharyngeal swab, or 1 ml whole blood in EDTA (purple top) tube, or tissue.

5. ADDITIONAL REQUIRED MATERIALS

- Disposable gloves
- Pipettes
- Sterile pipette tip
- Vortex mixer
- Centrifuge for microcentrifuge tubes
- Thermal cycler
- Electrophoresis kit
- UV transilluminator

6. PROCEDURE

Please read through the entire procedure before starting.

6.1 DNA PREPARATION

Various manufacturers offer DNA isolation kits. Please carry out the DNA isolation according to the manufacturers instructions. The following standard DNA Purification kit is recommended.

Product	Catalog No.	Manufacturer
Bioingentech™ Genomic DNA Purification Kit (50 test)	230040(50)	Bioingentech Biotechnology Inc.
Bioingentech™ Genomic DNA Purification Kit (100 test)	230040(100)	Bioingentech Biotechnology Inc.

6.2 AMPLIFICATION

1.- Prepare appropriate PCR Premix tubes and one PCR Premix tube for Positive control. Label.

2.- Add 6µl of DNase/RNase-free water into the PCR Premix tube to total volume as 11µl.

3.- Add 2µl of template DNA into the PCR Premix tube to total volume as 13µl.

4.- Add 6µl of DNase/RNase-free water and 2µl of PCR Positive control into a PCR Positive control Premix tube for monitoring of amplification and easy interpretation.

5.- Add mineral oil (11µl). This step is necessary, even when using a thermal cycler that employs a top heating method.

6.- Perform PCR reaction of samples as the below process, using a PCR thermal cycler.

PCR cycle		Temp.	Time
1 Cycle	Initial Denaturation	94°C	2 min.
	Denaturation	94°C	30 sec.
	Annealing	55°C	30 sec.
30 Cycles	Extension	72°C	30 sec.
	Final extension	72°C	5 min.

6.3 DETECTION OF AMPLIFIED PRODUCTS

1.- Prepare 1.5% agarose gel containing Ethidium bromide (Et-Br).

2.- Load 7µl of PCR product, 7µl of positive control and 2µl of Brig™ Molecular Weight marker on agarose gel without adding a loading-dye buffer and perform electrophoresis.

3.- Run electrophoresis by 100V (required about 30~40 minutes).

4.- Identify the result on ultra-violet (UV) transilluminator.

6.4 INTERPRETATION

- Expected PCR product size : 290 bp

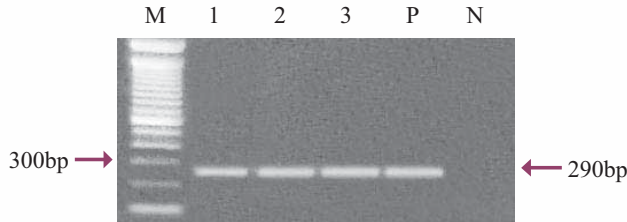


Fig1. Electrophoresis of PCR product by VetPCR™ EHV-1 Detection Kit
Lane M : Brig™ Molecular Weight Marker (Bioingentech Ltd.)
Lane 1~3 : EHV-1 Positive sample
Lane P : Positive control
Lane N : Negative control

7. NOTICE

- For research purpose only. Not for use in diagnostic procedures for clinical purposes. *For in Vitro Use Only.*
- Take care in handling of specimen to minimize risk of infection.
- The PCR process is covered by patents issued and applicable in certain countries. Bioingentech Biotechnology Inc. does not encourage or support the unauthorized or unlicensed use of the PCR process. Use of this product is recommended for persons that either have a license to perform PCR or are not required to obtain a license.

8. TROUBLE SHOOTING

- 1.- In the case of difficult to interpret results due to non-specific bands; reduce amount of template by 1/10 dilution, heated at 65° C for 5 min. and reacts again.
- 2.- Preparation of PCR reaction at room temperature may cause the non-specific band.
- 3.- All procedure should be carried out on ice.

9. ORDERING INFORMATION

Product	Catalog No.
VetPCR™ EHV-1 Detection Kit 48	VET0004E(48)
VetPCR™ EHV-1 Detection Kit 96	VET0004E(96)
Brig™ Molecular Weight Marker	24012



Bioingentech Ltd.

Salas 350, piso 2, Concepción, Chile
Telephone (56)-(41)-2790435
Fax (56)-(41)-2790435
info@bioingentech.com
www.bioingentech.com

Product use limitations warranty disclaimer

Bioingentech manufactures product for a number of intended uses. Please refer to the product label for the intended use statements for specific product. Bioingentech products contain chemicals which may be harmful if misused. Due care should be exercised with all Bioingentech products to prevent direct human contact. Each Bioingentech product is shipped with documentation stating specifications and other technical information. Bioingentech products are warranted to meet or exceed the stated specifications. Bioingentech's sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event products fail to perform as warranted. Bioingentech make no other warranty of any whatsoever, and SPECIFICALLY DISCLAIMS AND EXCLUDES ALL OTHER WARRANTIES OF ANY KIND OR NATURE WHATSOEVER, INCLUDING, WITHOUT LIMITATION, AS TO THE SUITABILITY, PRODUCTIVITY, DURABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, MERCHANTABILITY, CONDITION, OR ANY OTHER MATTE WITH RESPECT TO BIOINGENTECH PRODUCTS. In no event shall Bioingentech be liable for claims for any other damages, whether direct, incidental, foreseeable, consequential, or special (including but not limited to loss of use, revenue or profit), whether based upon warranty, contract, tort (including negligence) or strict liability arising in connection with the sale or the failure of Bioingentech products to perform in accordance with the stated specifications.

2004 – 2009 Bioingentech corporation. All rights reserved.

Products may be covered by pending or issued patents or may have certain limitations. Please visit our web site for more information.

All prices and specifications are subject to change without prior notice

Product claims are subject to change. Please contact Bioingentech technical services or access the Bioingentech online catalog for the most up-to-date information on Bioingentech products.